

EuroClone®

s@femate
EZ

Microbiological
Safety
Cabinet



EuroClone®
serving science through innovation

Stay safe the EZ way



Including premium features in an entry level priced cabinet, the new SafeMate EZ series allow everybody to stay safe in an easy way.

The "V"-shaped antiobstruction grill allows working without worrying about armrests, while the self positioning electrically operated front sash makes sure that the front aperture is always at the right size for comfort and safety.

Your Safety is our Commitment

No compromise for Operator, Product and Environment.

Protection guaranteed as required by EN12469:2000 standard.

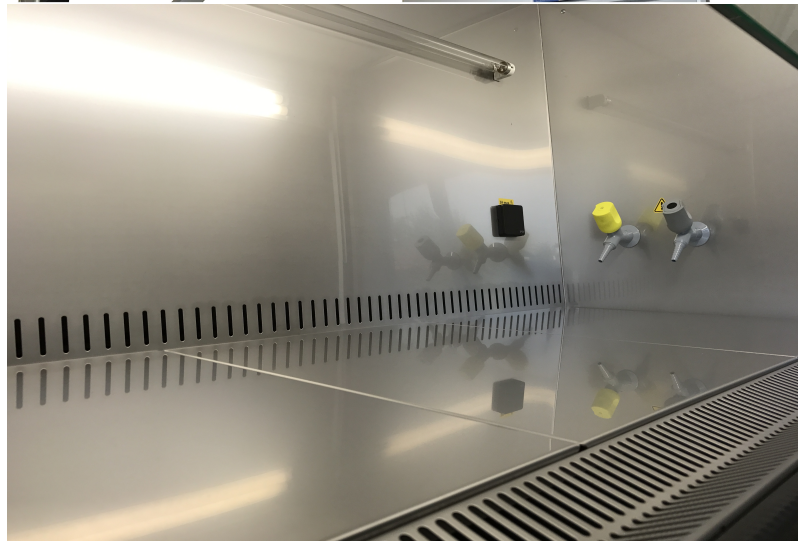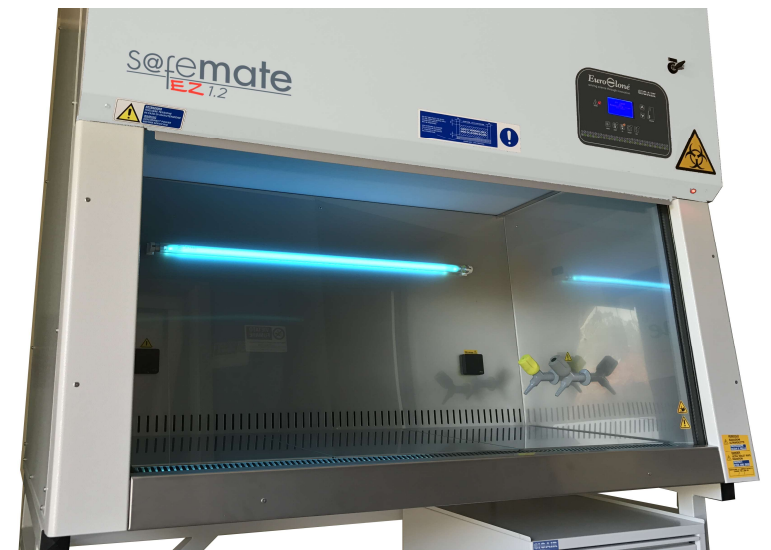



Feature rich!

Electrical front sash: the front glass is operated using the switches on the main control panel allowing effortless opening and closing of the working area.

V-shaped front grill: forget about armrests limiting your working position: SafeMate's V-shaped front grill ensures the front barrier is always at its best.

Customizable utilities: want more space in the working area? Do not take the optional taps if. Changed your mind and want the taps? Just buy the option and they will be installed in your cabinet even after sales!

Fully VHP compliant: with the optional VHP connector kit you can easily use any Hydrogen Peroxide vapour generation system to fully sterilize your cabinet.

Italian Design & Technology


Real Time
visualization of
airflows speeds

UV Lamp on back
wall (standard)


Electrically
operated front sash

Gas & Vacuum
taps (optional,
retrofitable)

Height adjustable
support stand

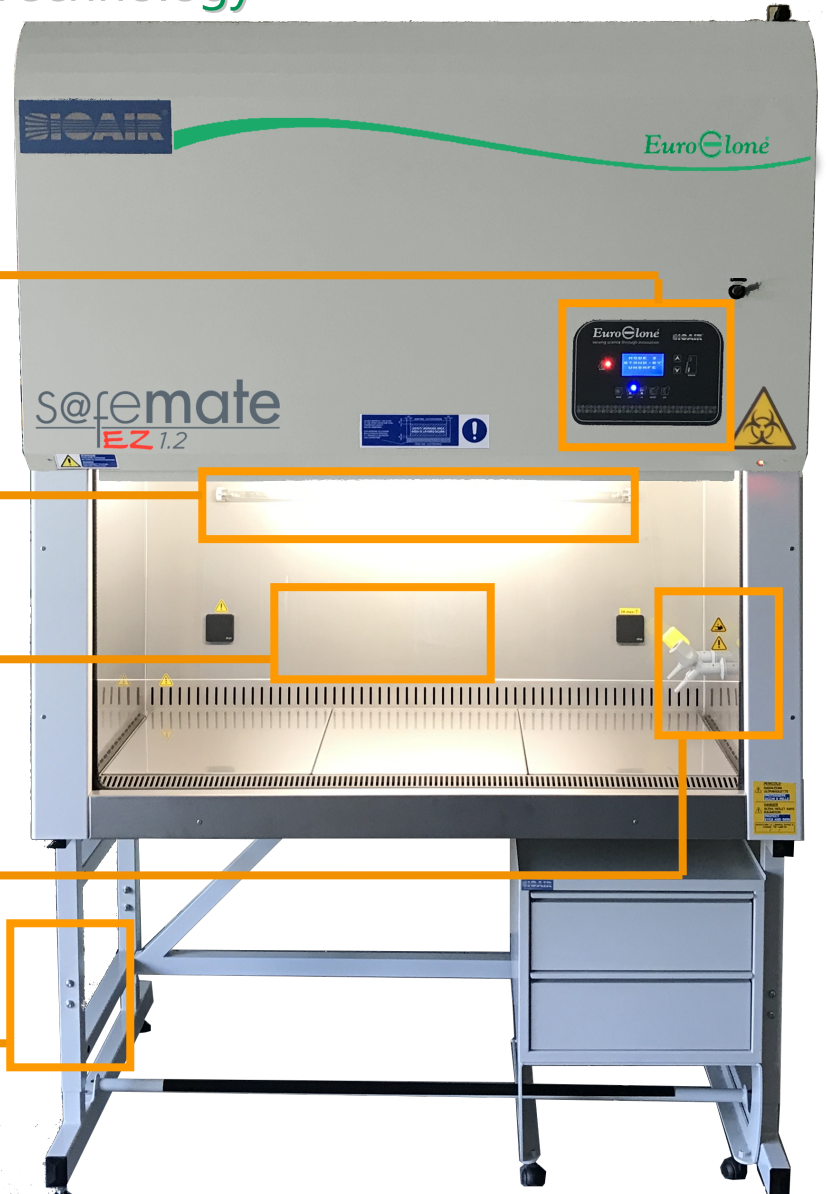

Choose your color!

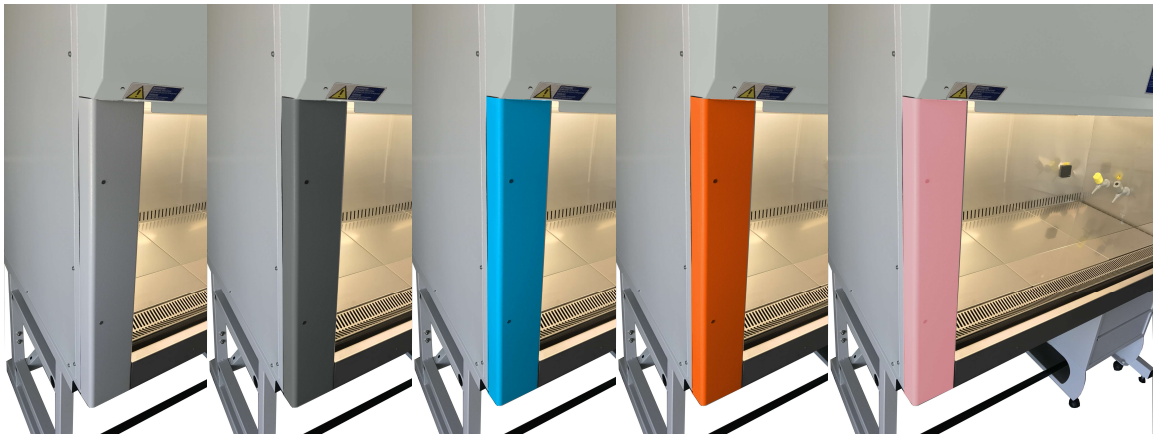

Our cabinets are completely made in Italy using components of Italian or European origins. We use only the best for our cabinets.

Main Specifications

State of the art AC motorblower enhances energy efficiency, reducing operating costs.

- Fully compliant with the EN 12469 safety standard
- Air and aerosol tight electrical sliding sash
- UV light on back wall
- 5° Sloping front aperture to maximise comfort.
- V-Shaped anti-obstruction front grill
- Optional utilities for gas and vacuum, installable in-field
- Real Time reporting of air speeds (inflow & downflow)
- Fully stainless steel working chamber and surface



Standard utilities

STANDARD ELECTRICAL EQUIPMENT	
Automatic electronic airflow velocity control PCB	✓
Motorblower (fan)	✓
Inverter	✓
Fluorescent lamps	✓
Sliding window electric motor	✓
Combustible gas solenoid valve.	✓

STANDARD UTILITIES	
Tap for combustible gas line	optional
Tap for inert fluids/vacuum line	optional
Auxiliary electrical service socket	✓
2 nd auxiliary electrical service socket	✓
UVC lamp socket	✓
Voltage-free contact (VFC) outlet	✓

Options & Accessories

PART No.	DESCRIPTION	NOTES	
AC10000	CHEST DRAWER	2 drawers – with castors	✓
AS1L410	SUPPORT STAND 1.2	Variable height	✓
AS1L610	SUPPORT STAND 1.8	Variable height	✓
AZ5K401	Combustible gas tap	Includes solenoid safety valve. Max 1 per unit	✓
AZ5K412	Vacuum/Air tap	Max 2 per unit (if no gas tap installed)	✓
AP6H100	Additional socket	Standard Schuko/Europlug	✓

Technical Data

DESCRIPTION	SIZE 1.2	SIZE 1.8
Part No. (solid work surface)	LDK4000	LDK6000
Part No. (perforated work surface)	LDK4001	LDK6001
SPECIFICATIONS		
Reference Standards:	IEC 61010-1:2010 / EN 61010-1:2010 IEC 61326-1:2012 / EN 61236-1:2013 EN 12469:2000	
Electrical insulating/protection class [IEC 61140]:	I	
Mains supply voltage:	220-240 V~ 50/60 Hz	
Required power line (W): (700 W service socket included)	1200	1750
Absorbed power (W): ^(*) (fan and light on only)	465	774
Window glass UVC radiations retention (%):	98	
Combustible gas fixture max pressure (mbar):	20	
Inert fluids/vacuum fixture max pressure (bar):	4	
Electrical service socket max current (A):	3	
WEIGHT AND SIZE		
Weight (kg):	256	360
Overall size L x D x H (mm): (without support stand)	1380 x 795 x 1450	1990 x 795 x 1450
Front aperture size L x H (mm):	1230 x 210	1840 x 210
Working space size L x D x H (mm):	1230 x 600 x 700	1840 x 600 x 700
MATERIALS		
Main structure:	cold rolled steel, stove enamel coated RAL 7035 + ABS plastics	
Working space surface:	stainless steel AISI 304 - 2B finishing	
Front window:	laminated safety glass	
PERFORMANCES		
Laminar Air Flow mean velocity [EN 12469](m/s):	0,33 ÷ 0,40	
Inflow Air Barrier mean velocity [EN 12469](m/s):	0,53 ±10%	
Exhaust Air flow rate (m3/h):	480 ±10%	600 ±10%
Exhaust Air flow ratio (%):	30 ±10	
Apf - Aperture Protection Factor [EN 12469]: (Retention efficiency at front aperture)	≥1,0 x 10 ⁵	
Working space air cleanliness class [EN 14644-1]:	ISO 5	
Illuminance [EN 12469] (lux):	>850	
Sound level [EN ISO 3744] (dB[A]): ^(**)	<56	<60
Vibration [EN 12469] (mm RMS):	<0,005	
Max increase inside cabinet in temperature from the ambient [EN 12469] (°C):	<5	
FILTERS		
Filters efficiency class [EN 1822-1]:	H14 ^(***)	
Filters global MPPS efficiency [EN 1822-1](%):	99,995	
MPPS diameter [EN1822-1](µm):	0,1 ÷ 0,3	

* Motorblower on, lights on (flow 0.28m/s, LED lights)

** Measured in operating conditions. Actual values at customer site may be different due to room structure

*** Efficiency higher than ULPA (Class F) as per IESP-RP-CC001

More than 40 years of experience

The experience of **EuroClone**[®] in manufacturing **Biohazard** and **Laminar Air Flow** cabinets goes back to the early 70s', when the brand *Gelaire*[®] became the "*gold standard*" for airborne contamination control in many laboratories throughout the world.

A family of **Recirculating Fume Hoods**, based on the adsorption of toxic vapors by means of charcoal filters, was successfully introduced a few years later, thus characterizing the Company as the only one really focused on the protection of the operators and inspired by its motto.

This unique know-how was cherished and brought to an even higher level of quality twenty-five years later, when under the name of **BioAir**[®], the entire range was completely re-designed to meet the growing requirements of the laboratory staff and the most stringent regulations.

At the top of the range, particularly noteworthy are the **Biohazard** (or Microbiological Safety) **Cabinets**, representing the sum of the Company's know-how certified to European standards (EN12469:2000) and complying with the Australian regulations; in other words, they are designed to provide the technicians with the maximum level of safety, when they are used according to GLP/GMP in their respective environments.

Today, in a plant occupying more than 2.800 square meters, **EuroClone**[®] manufactures a *complete range of microbiological safety cabinets, laminar flow cabinets and fume cupboards*, encompassing more than 15 models, with many of them available in different sizes; customized models and/or designed for specific applications can be produced thanks to the competence of a team of skilled engineers and dedicated workers.

The experience deriving from decades of sales and support to Cell Biologists, allowed EuroClone[®] to bring into the market an *extremely innovative CO₂ Incubator*, the **S@fegrow 188**, which is the result of a deep knowledge of the best conditions required by the most critical tissue culture methods, supported by the suggestions received from the scientists involved in growing cells *in vitro*.

The core business of the recently established **BioAir**[®] **Industrial** Team is the design, manufacturing and validation of customized equipment for the protection of the operator and of the product within *pharmaceutical and healthcare production facilities*.

This dedicated team will take advantage of the long experience and the production capacity acquired through laboratory LAF applications, to offer dedicated and complex equipment, ranging from **dispensing/sampling Downflow Booths** and **Clean Rooms**, to **RABS** and **Isolators** for Regenerative Medicine and Advanced Cell Therapy.





EuroClone®
serving science through innovation



EuroClone S.p.A.

Via Figino, 20/22 - 20016 Pero (MI) Italy
☎ +39 02 38195.1 - 📠 +39 02 38101465
✉ info@euroclone.it - www.euroclone.it

Quality Management Systems certified according to ISO 9001 and ISO 13485 international standards