

## Labo Gas Regulator

Quality – Made in Germany



For two different adjustable gas outlet pressures out of one gas supply


### Technical data:

- Outer dimensions:  
191 x 88 x 208 mm (w x d x h)
- Weight: 2.2 kg
- Gas inlet pressure: max. 3.0 bar
- Gas outlet pressure: 0 – 2.5 bar
- Use with CO<sub>2</sub>, N<sub>2</sub>
- Wall mounting possible (opening at the back)

REF 15971

### Characteristics:

- Necessary for the fluctuation-free operation of the Labo Gas Mixer (continuous gas consumption) if the same gas source is used to provide gas for CO<sub>2</sub> incubators (pulsed gas consumption)
- Adapted for gas supply by gas bottle with gas pressure reducer or facility installed gas supply
- Suitable for operation in various environments where different gas pressures need to be generated from one supply source

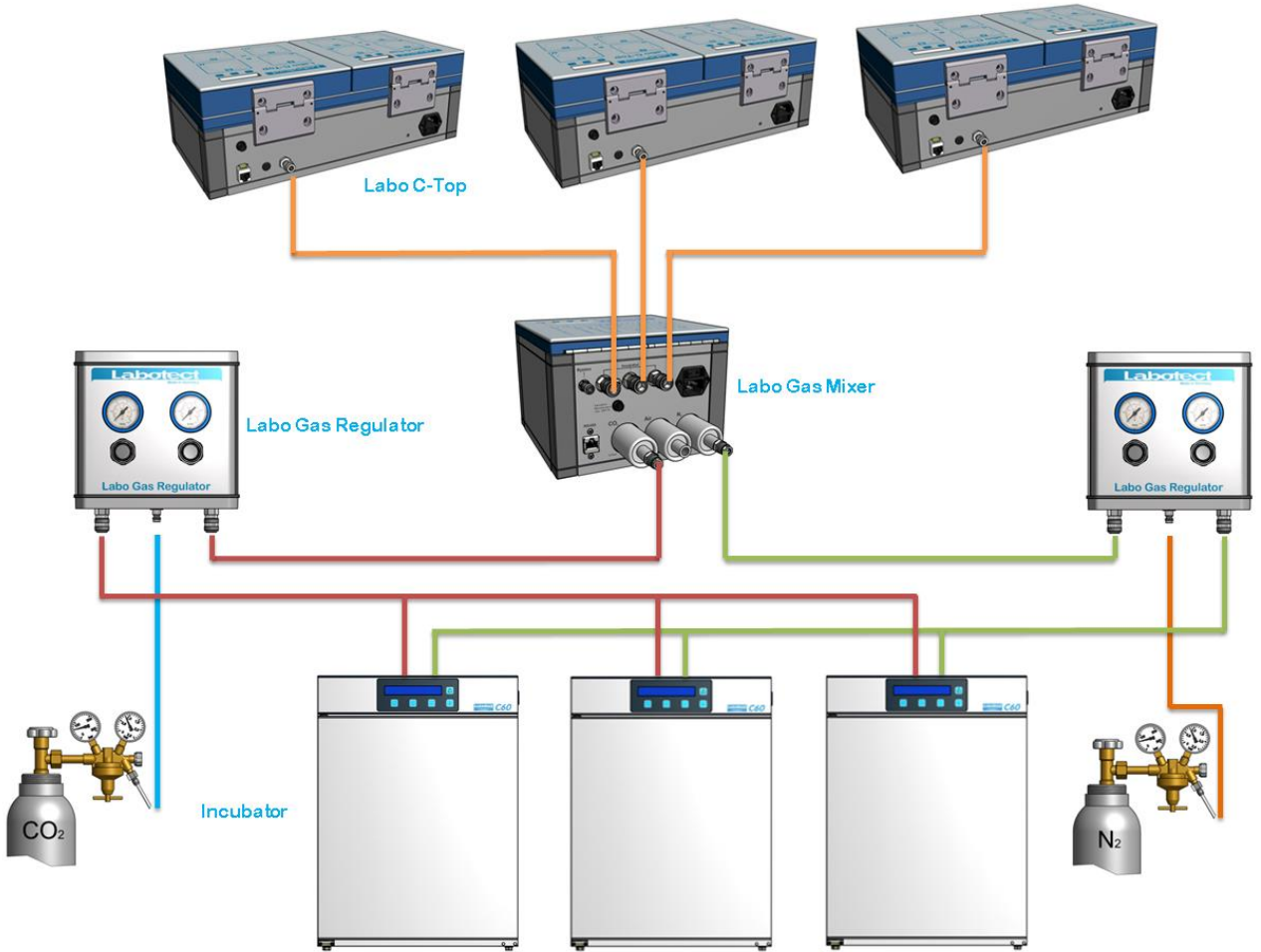

Picture: Application example of Labo Gas Regulator with Labo Gas Mixer, 3 x Labo C-Top and 3 x C60 incubator



Rev. 2\_01/2017