



EDICAL-PHARMACEUTICAL REFRIGERATOR MPR 1160 xPR



- Structure: upright type, realized in sheet steel epoxy coated white colour both in- and outside or in stainless steel 18/10 AISI 304 both in- and outside. All the internal corners are rounded to make easy any cleaning operations and the internal bottom is tray type to contain
- Insulation: high density (40 Kg/m³) foamed-in-place polyurethane, with a thickness of 75 mm. CFC-free
- Feet: Nr 4 made in stainless steel 18/10 AISI 304. adjustable in height for levelling
- Glass doors: Nr 2, hinged, side-by-side, made a perimetrical aluminium frame and three layers of tempered insulating glass, with vacuum among the glass layers in order to increase the insulating rate. The doors are provided with a removable magnetic rubber gasket in order to grant a perfect hermetic closing. The handles are realized in aluminium, integrated in the door profile, fullheight for an higher comfort of use. The doors are also fitted with a spring loaded automatic closing device for openings inferior then 90° and special door switches that stop the internal ventilation at the door opening (to limit the air cold loss)
- Internal equipment: Nr 6 (3 for each door) storage open wire shelves realized in sheet steel with a strong plastic coating (rust-proof material). The shelves are mounted on special anti-tilt stainless steel slides that allows the shelf extraction and they can be easily removed and adjusted in height without the use of any tool. The special internal racking-system allows a great flexibility of the internal fitting with the possibility of interchange and combine both

shelves and drawers (drawers optionally available)

- Shelf dimensions (W x D cm): 53 x 55
- Shelf loading capability (Kg): 35 (with uniformly stored material)
- Internal lighting: Nr 2 LED tubes, mounted in the side walls, with activation through a special switch located on the control panel. They grant energy saving and prevent the internal heating produced by the traditional neon tubes
- Control panel "xPRO": located in the top part of the structure (above the door), microprocessor operating, with capacitive touch-buttons and colour display. The controller supervises all the appliance functions through the EEN (EVERmed Easy Navigation) interface, with icons and written messages, allowing a very intuitive interactivity between the appliance and the user, for an easy setting of all appliance functionalities according to needs and the control of the workings status. The control panel is battery operating in order to grant the functioning of the display and alarms even in case of power failure. The main functions panel are:
 - High resolution 3,5" TFT colour display, LED backlit, showing the appliance working status (actual temperature, compressor functioning, defrosting, alarm running, etc.) and allowing the surfing through the various controller menu for setting and checking purposes
 - Digital temperature displaying and adjusting with an accuracy of 0,1°C
 - Keyboard touch-buttons with automatic block safety protection, to avoid tampering from unauthorized persons
 - Dual-level safety access, with PIN code protections, on both ON/OFF switch and keyboard touch-buttons. Reprogrammable PIN codes for increased security from access to the unauthorized personnel
 - Back-up battery, with automatic recharging, granting till 48 hours autonomy for the control panel functioning in case of power failure. The charging status is continuously shown on the display for immediate checking
 - Visual and acoustic alarm signalling (with automatic resetting) for:
 - high and low temperature, with limits programmable from the user through a dedicated menu
 - door ajar, delayed of 2 minutes to allow the standard operations. Programmable alarm delay through a dedicated menu









EDICAL-PHARMACEUTICAL REFRIGERATOR MPR 1160 xPR

- power failure
- battery discharged
- condenser dirty
- anti-freezing evaporator
- sensors failure
- Muting facilities for the acoustic alarms, with continuous warning indication on the display and automatic ring-back if the alarm situation persists. Programmable ring-back delay through a dedicated menu
- Test facilities for both control panel functioning, battery charging status and alarm limits
- Alarms log history for the latest 30 alarm conditions (not deletable), with the possibility of verify on the display, the alarm condition and all the related informations (starting time, ending time, highest/lowest peak of temperature reached where applicable)
- Door openings history for the latest 30 days (not deletable), with the possibility of verify on the display, the number of door openings, how many giving alarms and the total time that the door is kept open day-by-day
- Defrosts history for the latest 30 days (not deletable), with the possibility of verify on the type of defrost, starting time and the duration
- Internal ventilation setting facilities, in order to optimize the air-circulation according to the appliance using conditions and needs
- Real-time temperature chart, with surfing possibilities of the previous 52 hours of operations and zooming facilities for a better control of a period of interest
- Written indications on the display selectable from 3 languages (Italian, English and French)
- Pt1000 sensors for a better accuracy of the temperature control
- NO/NC contact to remote the alarm signals
- Software upgrading facilities
- Cooling unit: top mounted, with the condensing unit compounded by Nr 1 hermetic compressors and Nr 1 finned condenser, air cooled through a fan. Also, in the ceiling of the storage chamber, is mounted the finned evaporator with copper pipes, that is housed into a stainless steel cover that incorporate the fan. All the mounted components are industrial grade to grant the maximum reliability
- Refrigerant: R404a CFC-free
- Refrigeration: forced-air, through a fan, granting the maximum temperature uniformity and stability inside of the
- Defrosting: completely automatic, thermostat controlled. The condensate water is automatically channelled into an heated tray, located under cabinet, for the automatic evaporation of the condensate water
- Temperature range: infinitely adjustable between +2°C / +15°C
- Voltage (V/ph/Hz): 220-230/1/50
- Plug: Schuko type
- Breakers: Nr 2 glass cartridge fuses with a rating of 10A, at protection of the appliance
- Noise level (dB(A)): ≤ 52 Gross capacity (litres): 1295 Net capacity (litres): 1160
- **Dimensions (W x D x H cm)**: 150 x 70 x 210
- Net weight (Kg): 230
- Packed dimensions (W x D x H cm): 158 x 98 x 224 (3,47 m³)
- Gross weight (Kg): 260









EDICAL-PHARMACEUTICAL REFRIGERATOR MPR 1160 xPR

ACCESSORIES AVAILABLE



Extra wire shelf



Wire basket



Plastic partitions for drawers



Kit of castors



Temperature chart recorder



Data logger function with USB port



Internal electrical socket



Internal spark-proofing




GSM/GPRS telephone dialer




Wooden crate packing



Extra stainless steel wire shelf

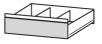

Stainless steel drawer

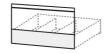

Front top for drawers



Safety door lock with key



Electronic temperature printer




Ballasted sensor for simulation of product temperature alarm



Cable port with cover



Back-up cooling unit

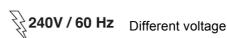



^{*} The declared performances are referred at +35°C ambient temperature, relative humidity of 70% and without any internal thermal mass (empty

^{*} Should it be necessary or useful without undermining any model essential features, technical and stylistic characteristics are subject to be changed without manufacturer's previous notice.



