

ahn myLab VT-02

Vortex Mixer 2800 rpm

Usage - everyday mixing operations

Maximum Speed - adjustable up to 2800 rpm

Comfortable use - selection of different modes (touch, stand-by and continuous), depending on the desired use / smooth and silent operation


Stable operation - through inbuilt counter balance system

Multiple applications - optional attachments for micro plate, micro tubes and flask (*only applicable in combination with the universal platform 7-110-01-0) /sticker for permissible speed for various attachments

Small footprint saves valuable bench space - suitable for use in hoods in the temperature range from 5 to 40°C



Optional attachments*

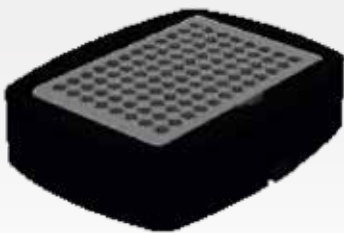

Foam for micro plate



Foam for micro tubes



Rubber holder for flask

Ordering information

Description	Cat. No.
AHN myLab® VT-02 Vortex Mixer 2800 rpm	7-110-00-0

Specifications	
Orbital diameter	4 mm
Speed accuracy	± 5.0%
Max. loading capacity	500 g
Protection class	IP21
Power supply	220 240 VAC, 50 Hz
Dimensions (WxDxH)	205x136x138.5 mm
Weight	3.0 kg

Delivery Package
1 pc. Vortex Mixer 2800 rpm
1 pc. Power supply adaptor
1 pc. Standard cup head (pre-installed)
1 pc. Instruction manual
1 pc. Warranty card

Optional attachments *	Cat No.
Universal platform	7-110-01-0
Foam for micro tubes*	7-110-02-0
Foam for micro plate*	7-110-03-0
Rubber holder for flask*	7-110-04-0

