

ahn myLab VT-03

Digital Vortex Mixer 4200 rpm

Usage - comfortable everyday mixing operations

Maximum Speed - adjustable up to 4200 rpm

Digital display - adjustable speed setting with an interval of 10 rpm and time setting from 0 - 999 minutes / displays speed and time alternatively at every 5 seconds

Stable operation - through inbuilt counter balance system

Robust design - brushless DC motor enables maintenance free and long lasting operation


Programmable pulse mode - ON time and OFF time can be individually modified

Multiple applications - optional attachments for micro plate, micro tubes and flask (*only applicable in combination with the universal platform 7-110-01-0) / sticker for permissible speed for various attachments

Small footprint saves valuable bench space - suitable for use in hoods and cold rooms



Optional attachments*

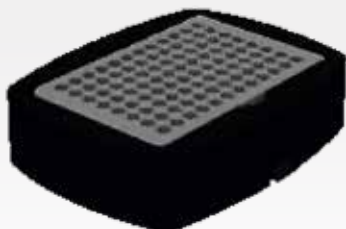

Foam for micro plate



Foam for micro tubes



Rubber holder for flask

Ordering information

Description	Cat. No.
AHN myLab® VT-03 Digital Vortex Mixer 4200 rpm	7-120-00-0

Specifications	
Orbital diameter	4 mm
Different modes	touch, continuous
Speed setting	300-4200 rpm
Max. loading capacity	500 g
Dimensions (WxDxH)	205x136x138.5 mm
Weight	2.8 kg

*only applicable in combination with the universal platform 7-110-01-0.

Delivery Package
1 pc. Digital Vortex Mixer 4200 rpm
1 pc. Power supply adaptor
1 pc. Standard cup head (pre-installed)
1 pc. Instruction manual
1 pc. Warranty card

Optional attachments *	Cat No.
Universal platform	7-110-01-0
Foam for micro tubes*	7-110-02-0
Foam for micro plate*	7-110-03-0
Rubber holder for flask*	7-110-04-0

